

# Water Baths



Speciality - Circulating & Deep Chamber



SHEL LAB special application water baths boast the same features as our other high-performance models, such as: being highly accurate, easy-to-use, rugged, and safe. SHEL LAB speciality water baths include a polymer coating for quick cleaning.

SHEL LAB was the first to introduce the non-contact recessed heating element design to the analytical research market. This curtails element burnout and eliminates tank hot spots that are chronic challenges for other water baths.

The SHEL LAB digital water baths remain an industry standard for precise control and quality of design. The operator is able to calibrate this bath using the touch pad controls. This is a convenient feature for facilities concerned with IQ/OQ/PQ validation. The SHEL LAB constant temperature digital water baths are perfect for conducting a host of applications including: bacteriological examinations, food processing/QC procedures and microbiology assays.


## FEATURES

- Easy-to-Clean Stainless Steel Tank
- Pocket Handles for Easy Lifting
- Recessed Heating Element Prevents "Burnout"
- Unique SHEL LAB Design Eliminates "Hot Spots"
- Thermometer Clip Included

*SHEL LAB digital water baths are durable and dependable.*

## Water Bath Applications

- Sample Thawing
- Bacteriological Examinations
- Warming Reagents
- Coliform Determinations
- Microbiological Assays

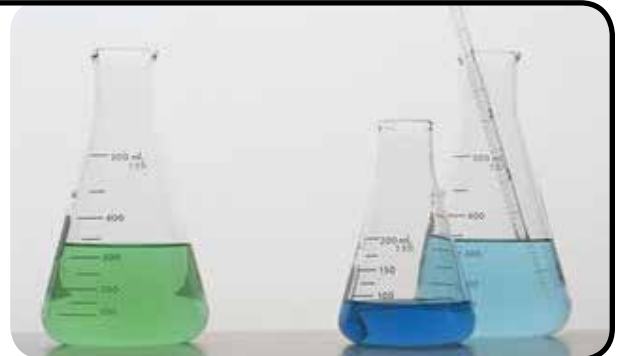

## Precise Temperature Control -

- Over Temperature Protection
- Temperature Uniformity +/- 0.2% at 37°C
- Temperature Range Ambient + 5°C to 80°C

**2 Year  
Limited  
Warranty!\***



Speciality Water Baths	Model Number		
	110-120V 220-230V Chambers	SWBC22 SWBC22-2 Circulating	SWB30 SWB30-2 Deep Chamber
Exterior Dimensions (wxdxh)	Inches	15.0 x 25.5 x 9.0	16.0 x 19.3 x 15.3
	cm	38.1 x 64.8 x 22.9	40.7 x 48.9 x 38.8
Bath Tank Dimensions (wxdxh)	Inches	11.5 x 19.5 x 6.2	12.0 x 14.0 x 12.0
	cm	29.2 x 49.5 x 15.8	30.4 x 35.5 x 30.4
Tank Capacity	L	23	30
Stainless Steel Gable Cover(s)	Yes/No	Yes - One Included	Yes - One Included

\*Tank is excluded from the 2 year warranty.

Specifications are based on nominal values at an ambient temperature of 25°C and a line voltage of 120V or 240V respectively. The temperature data is determined in accordance with DIN 12880. We reserve the right to alter technical specifications at any time.

For a complete list of product specifications, call (800)322-4897 or visit the SHELL LAB website, [www.shellab.com](http://www.shellab.com)

0740534 6/14



[www.shellab.com](http://www.shellab.com)  
1-800-322-4897